

sacred[®]
Symbol of quality



DENTAL, ORTHODONTIC &
IMPLANT SURGERY INSTRUMENTS

www.sacred.dental
info@sacred.com.pk



INTRODUCTION

Sacred International specialises in Dental, Orthodontic and Implant Instruments. With more than 35 years of experience, we have been manufacturing quality dental products which are aesthetically designed and easy to use. As manufacturers we produce high end instruments which are crafted with the utmost attention to details.

Our company maintains high levels of innovation, proficiency, technical prowess and constant quality checks to deliver integrated solutions that work all around the world since 1987 specially in Europe, USA, Japan, Turkey, Korea and many more countries around the globe. We are constantly developing new and innovative products and improving the current ones. Our customers feedback has helped us achieve in being more effective and efficient over time and Sacred International is now known as a Symbol of Quality and Reliability.

Thanks to our team members which consist of highly skilled, trained and efficient employees who work around the clock to help turn our visions into reality. With the reliability and quality perfected products we have developed loyal customers and hope to expand this family around the globe.

MAINTAINING HIGH LEVELS OF INNOVATION,
PROFICIENCY, TECHNICAL PROWESS
AND CONSTANT QUALITY CHECKS.



VISION

Our vision set out nearly 35 years ago was to become essential to our customers & advance in the dental industry by providing efficient, innovative designs and a highly reliable customer service which would not only make it easier for our customers to communicate but helps us to develop a bond of trust – something that we never compromised on. Now we proudly stand by and continue working on this ideology set years ago alongside tapping in the modern era through revolution and modern technology. Creating an environment where we deliver solutions in the shape of ergonomically designed dental instruments with respect to Dental Professionals requirement.



MISSION

A perfect combination of traditional and modern methods is used by our skilled technicians to manufacture high-end quality. Our mission highly relies on customer satisfaction as well. Therefore, highly trained representatives work around the clock for customer care service who are ready to assist you with any problem or query. We believe in building and strengthening our bond with our customers and continuously expand this family based on quality work and prompt deliveries without compromising the trust of valued customers.


CERTIFICATIONS

Sacred International operates under compliance with exact requirements of Directive 93/42/EEC, is certified by:

- FDA (USA)
- JPAL (JAPAN)
- KFDA (KOREA)
- MHRA (United Kingdom)
- EUDAMED (European Union)

And our quality management system is ISO certified in accordance with:

- ISO: 9001:2015
- ISO: 13485:2016

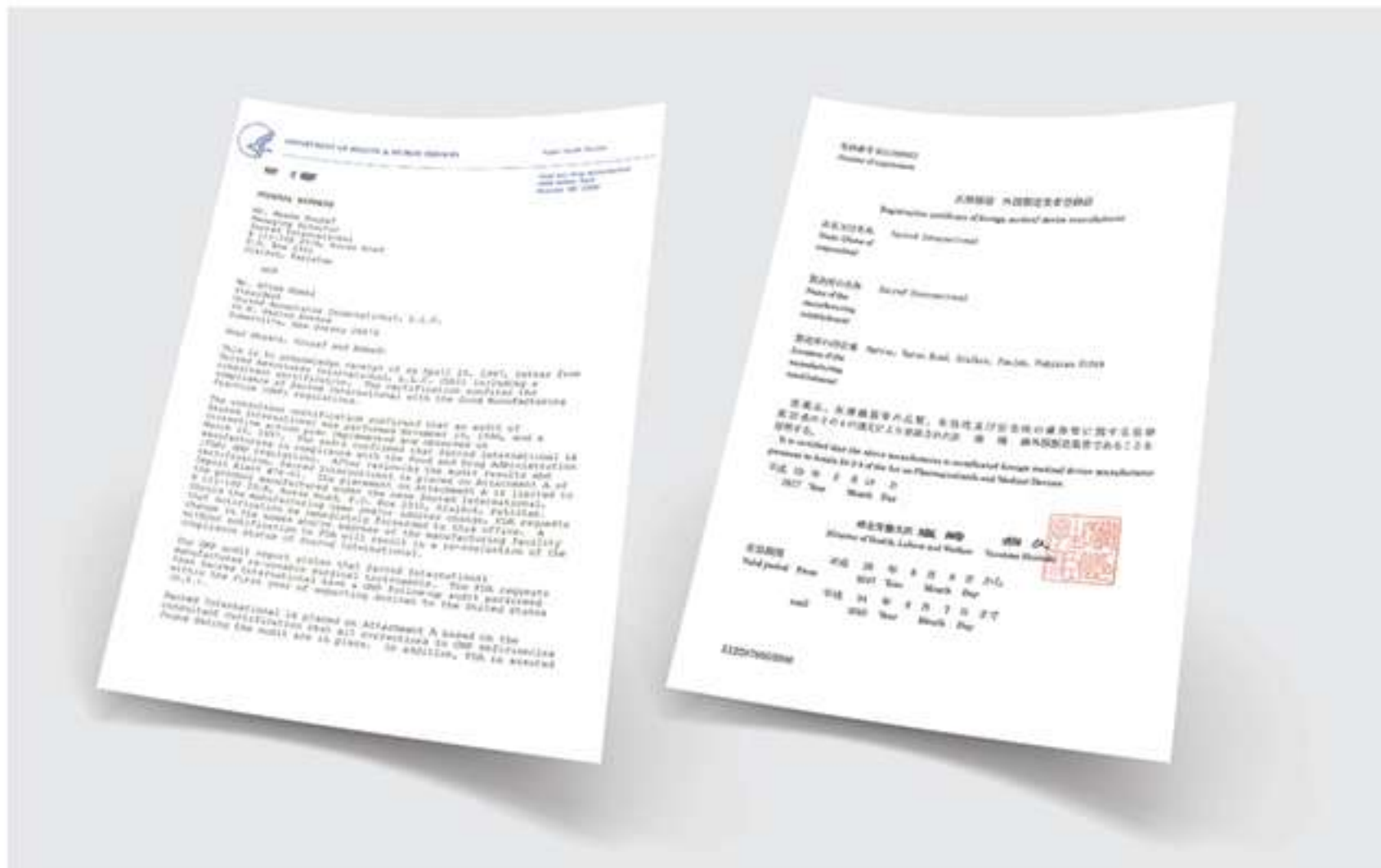




CLEAR ALIGNER PLIERS



OVERVIEW:

The Aligner Pliers help in recovering the teeth that may have been affected whilst an ongoing treatment or stubborn teeth that haven't been affected during the treatment.

The Aligner Pliers help in creating smaller indents which helps in the movement of the teeth by applying force and it helps in increasing retention. These aligners are much easier and sustainable in eliminating prolonged treatments and helps in torque movements, rotation, lingual movements for an increased retention on teeth.

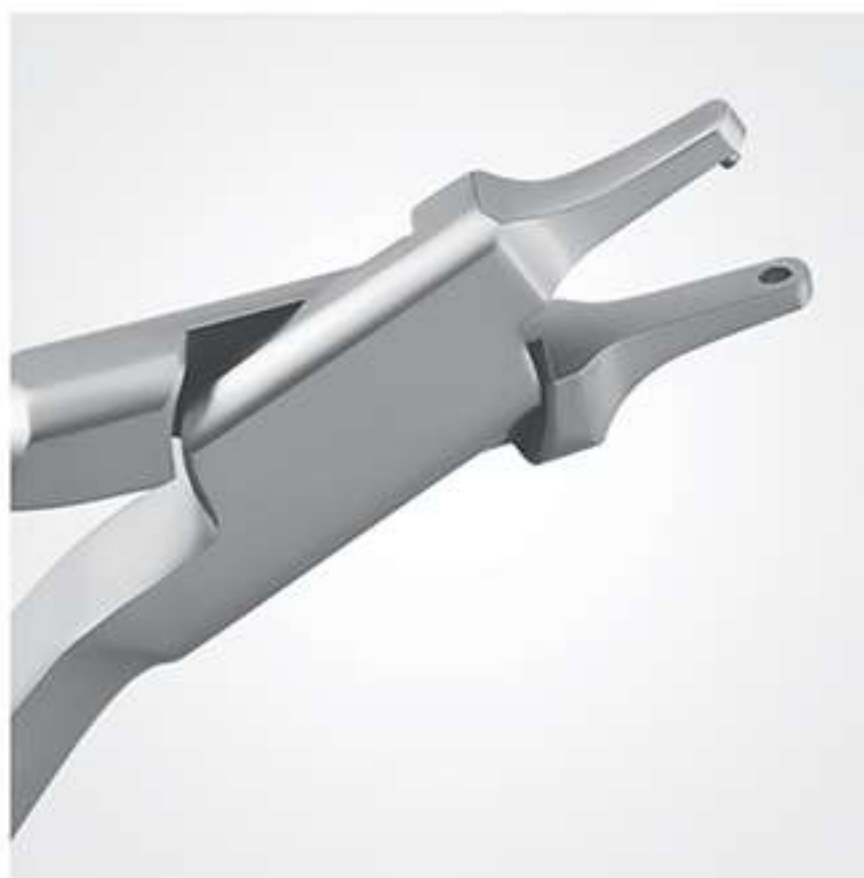

SPOT PLIER

AP001

The Spot Plier helps by creating small, rounded depressions for a precise contact which help to accent in creating movement of small individual teeth and for adding retention. This is normally used in underlining the movements of any tooth placed in the aligners that need a slight push.





HOLE PUNCH PLIER

AP002

The Hole Punch Plier creates 6.0mm semi-circular cut outs for an ease for the clear aligners around bonded buttons used for elastics. It precisely also nip's away plastic to provide soft tissue clearance to prevent aligner impingement of tissue.



HORIZONTAL PLIER

AP007


The Horizontal Plier helps in creating an indent to accentuate individual root torques. It also amplifies the retention of retainers or clear aligners by placing an indent into tooth undercut bonded attachment.



TEAR DROP PLIER

AP008

Tear Drop Plier craves at the gingival margin of clear aligners for where elastic hooks are required which provides a reservoir to hold the elastics in place as the aligner is seated which permits an easier manipulation of elastics by the patient.

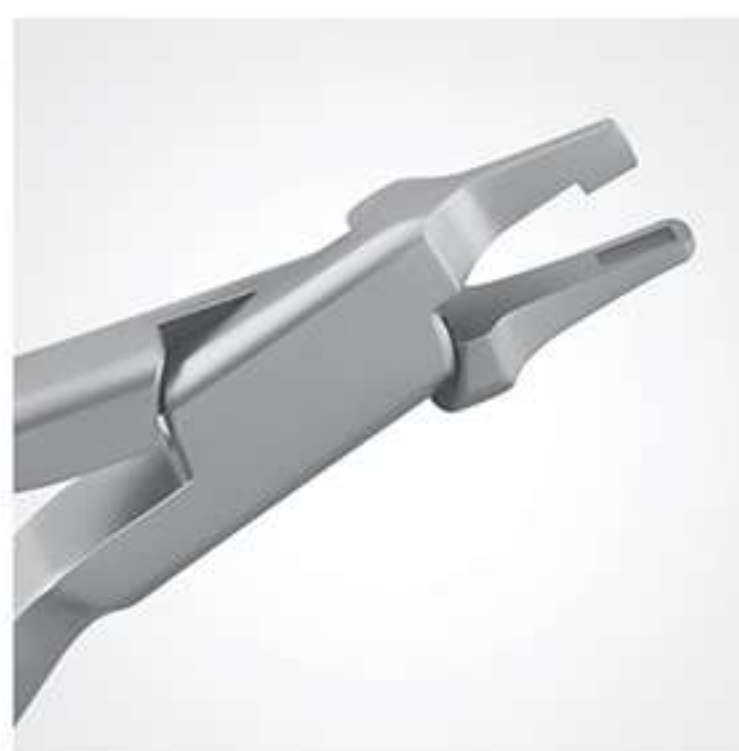

VERTICLE PLIER

AP009

Vertical Pliers help in rotational couple that can also be created using vertical indents on facial and lingual surfaces on parallel sides of a tooth. The Vertical plier can also be used to produce an indent to produce a mild cut to enhance resolution of in-and-out discrepancies.



POINT ELEVATION PLIER (1.5 mm)



AP006

1.5 mm Point Elevation Pliers play a vital role for creating punctiform elevations in the aligner. They allow a special activation which increases the pressure on the individual teeth.



FLEXIBLE PERIOTOMES

Flexible Periostomes is also one of our most enjoyed products among our customers. It's made up of high-quality stainless steel. It helps in cutting periodontal ligaments for atraumatic dental tooth extractions and dental implant placement. It has sharp thin blades which help in the extraction and removal of teeth with very less to no damage to the surrounding alveolar bone.

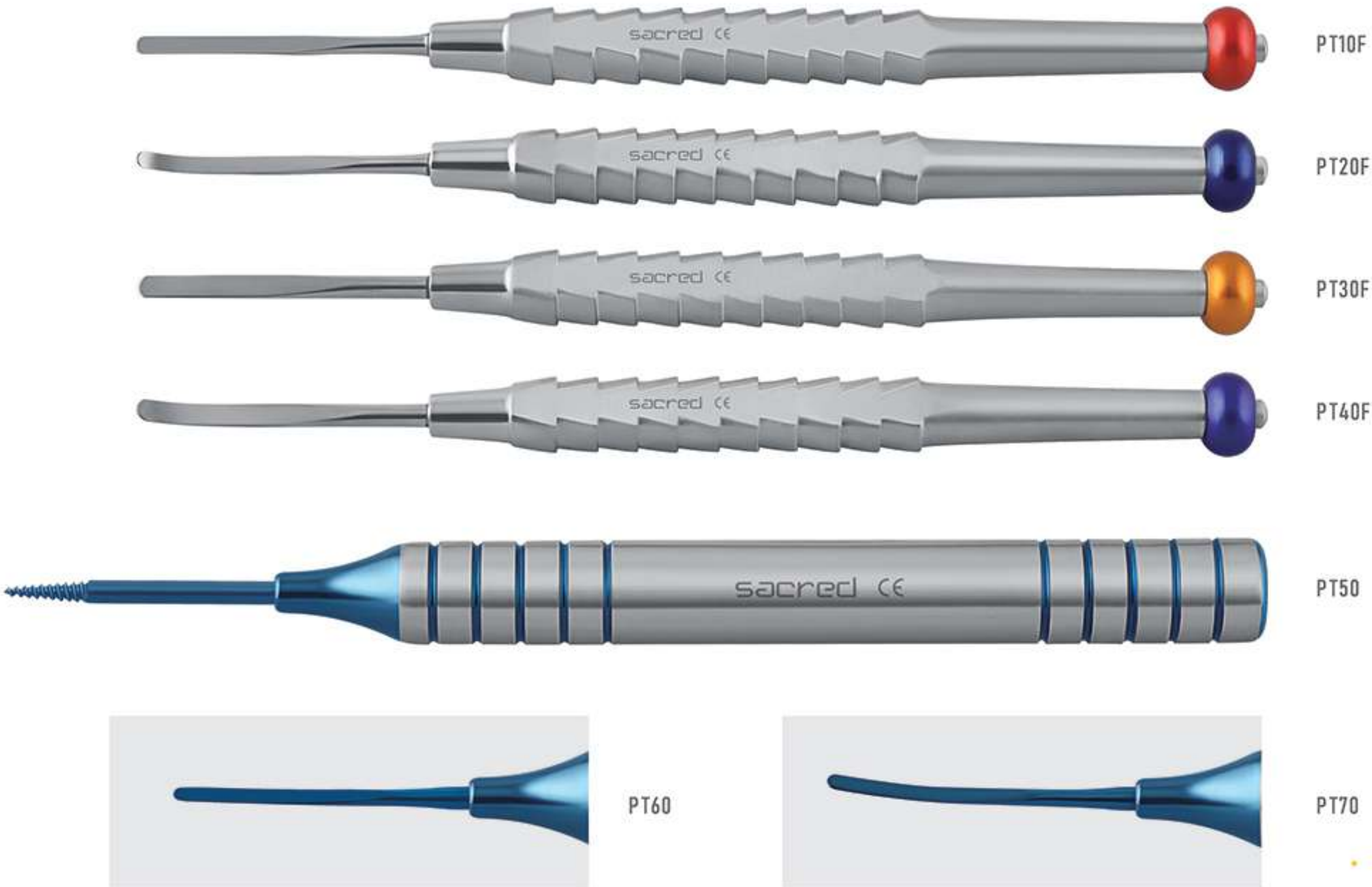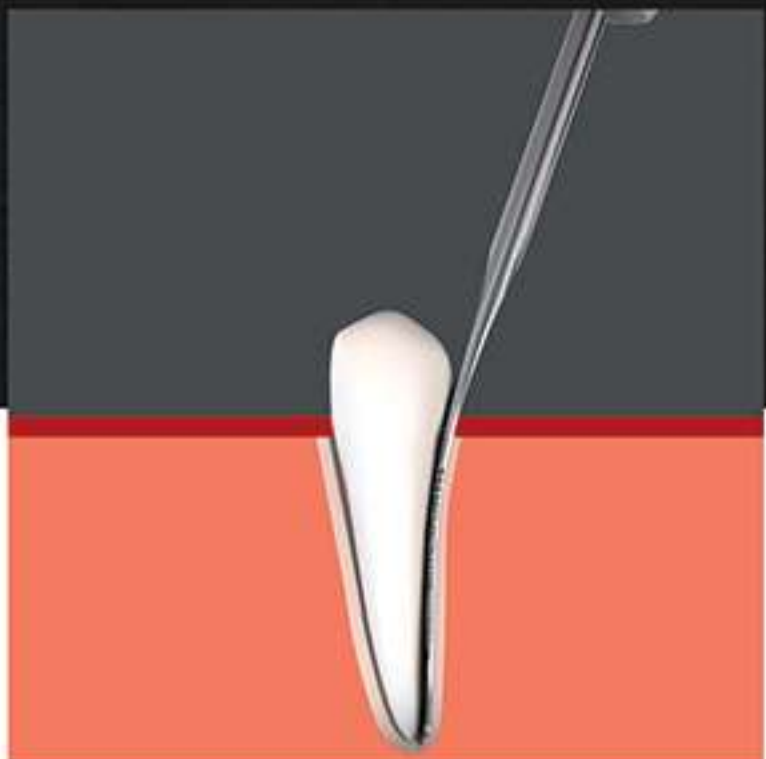

RUBBER DAM KIT

sacred
Symbol of quality



Introducing new Re-designed & Technically Altered Rubber Dam Kit that helps clinical operations to be carried out efficiently with ease.

FILLING INSTRUMENTS - PRO EDITION



025-074N

025-069N

025-083N

025-108N

025-123N

SCALPEL HANDLES WITH SILICONE GRIPS

The new Scalpel Handle comes with distinct features making itself a best choice for dental professionals.

- Reduced 8mm diameter.
- 45% decreased volume to ratio.
- More fun in crafting glowing smiles with a vast range of personalised colours.


www.sacred.dental

THE ACHROMATIC™ COLLECTION



The ACHROMATIC™ Collection is mapped out keeping in mind the clinical requirements to further ease the use of instruments during surgery.

- Distinct ceramic black coating with aesthetic matte visuals.
- Smooth performance.
- Reduced light reflection.
- Maximum comfort and control while performing surgery.
- Note: Also comes in black titanium finish to further enhance the durability and corrosion resistance.

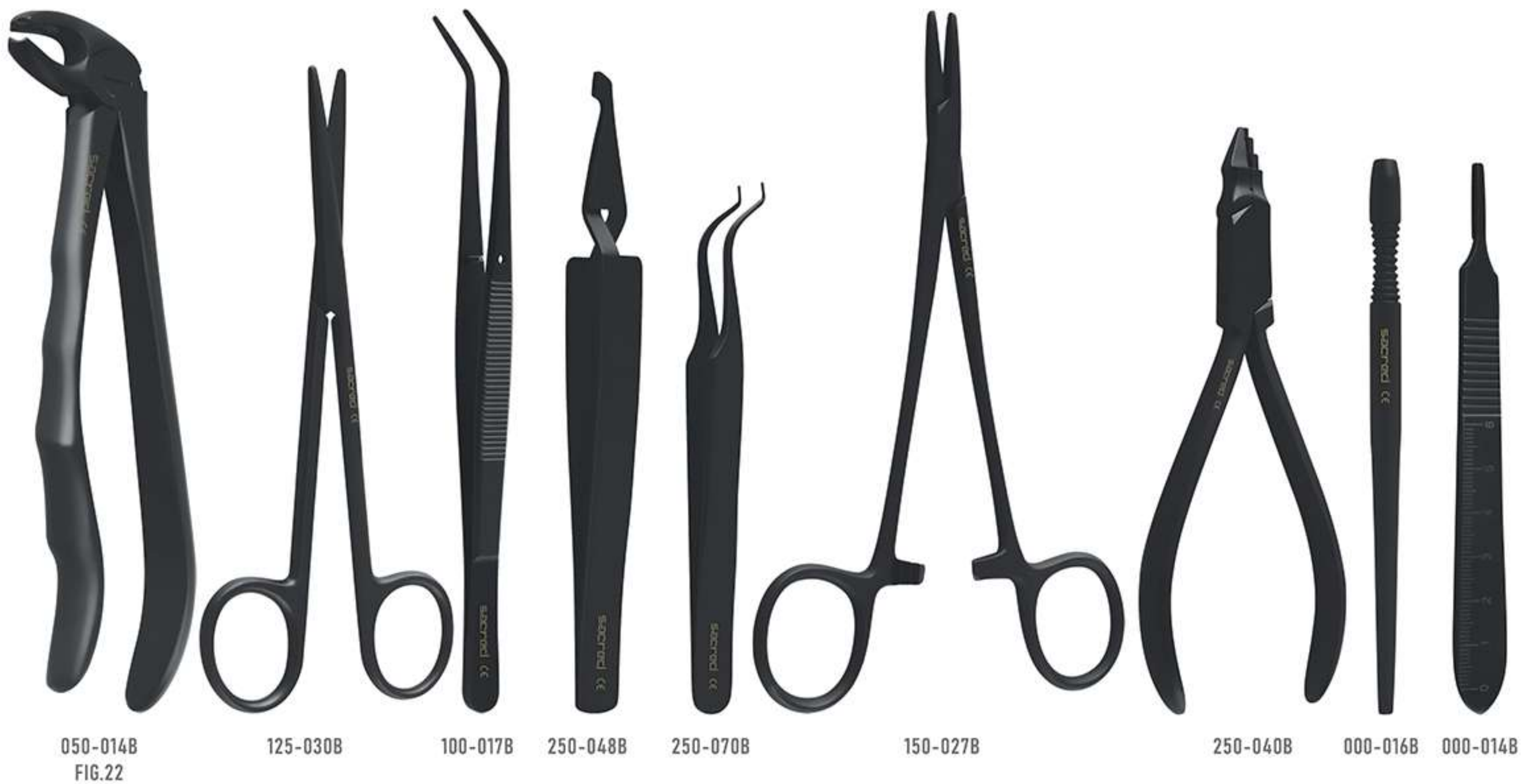

Ultra-thin, flexible composite instrument which acts as a solid brush by integrating Nickle Titanium material which is known for Its super elasticity properties.

DENTAL LUXATORS

sacred[®]
Symbol of quality




THE DENTAL LUXATOR TAKES UP DENTISTRY TO ANOTHER LEVEL.

- Sharp and less concave blade.
- Softer steel than elevator.
- Separates perforating fibres within the periodontal ligament.
- Slack the tooth prior to extraction.





ORTHODONTICS PLIERS

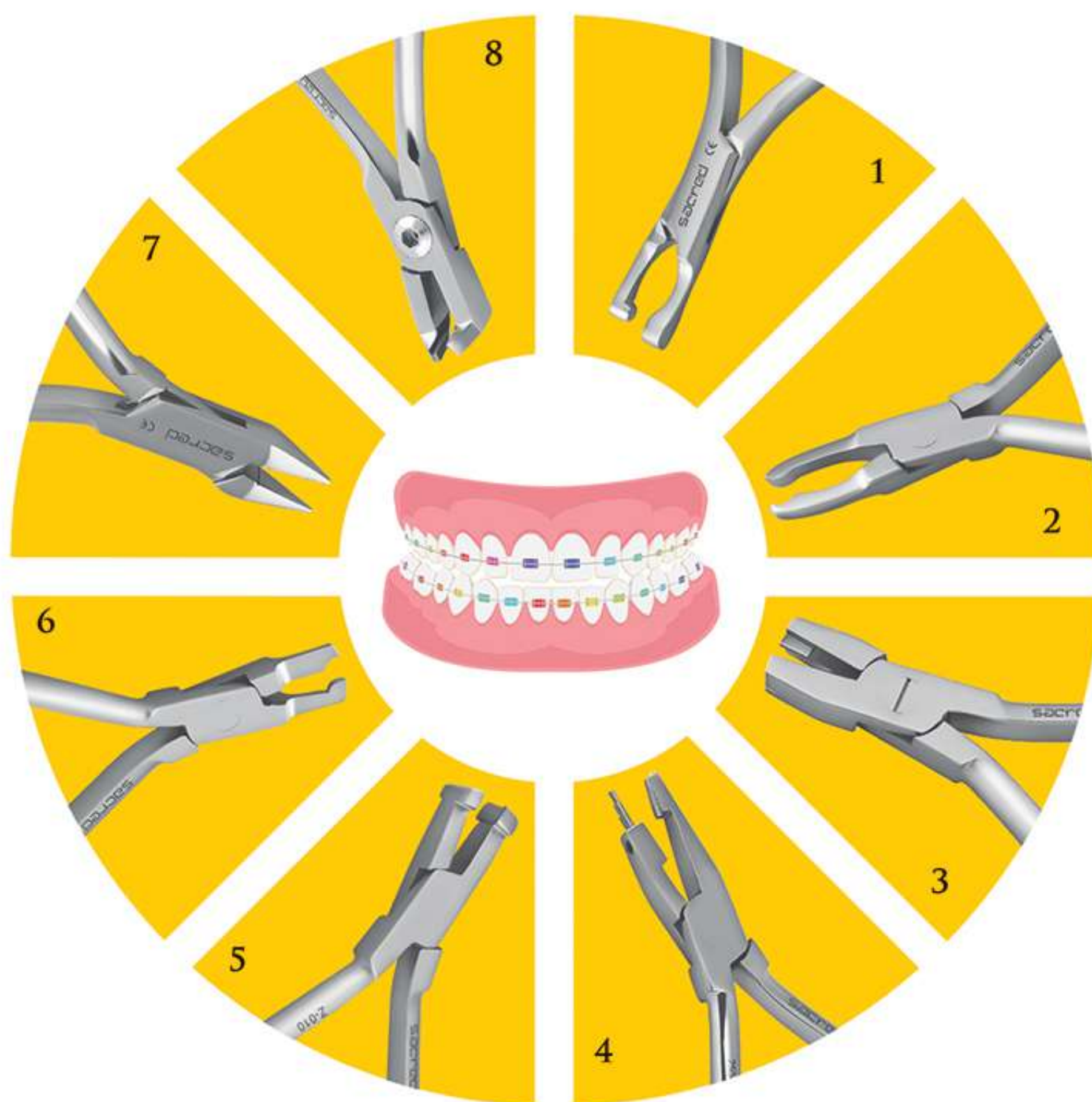

- | | |
|-----------------------------|------|
| 1. BALL & SOCKET PLIER | 060B |
| 2. JOHNSON BAND PLIER | 061B |
| 3. NANCE TIE BACK PLIER | 062 |
| 4. TWEED LOOP FORMING PLIER | 063B |
| 5. Z FORMING PLIER | 064B |
| 6. BAND CRIMPING PLIER | 065B |
| 7. THREE GROOVE PLIER | 066B |
| 8. DISTAL END CUTTER | 001 |

TORQUE BENDING PLIER



- THEY ARE USED IN PAIRS FOR APPLYING LABIAL OR LINGUAL TORQUE TO ARCHWIRES.
- PARALLEL GRIP PREVENTS CRUSHING OF WIRES



ACCURATE PLACEMENT OF TORQUE
TO ARCHWIRES


SCISSORS Premium

sacred
Symbol of quality



125-004P

Introducing an addition to our **ACHROMATIC™** Collection, Premium Scissors SuperCut with the help of tungsten carbide inserts in the blade & HD black coating on the surface.

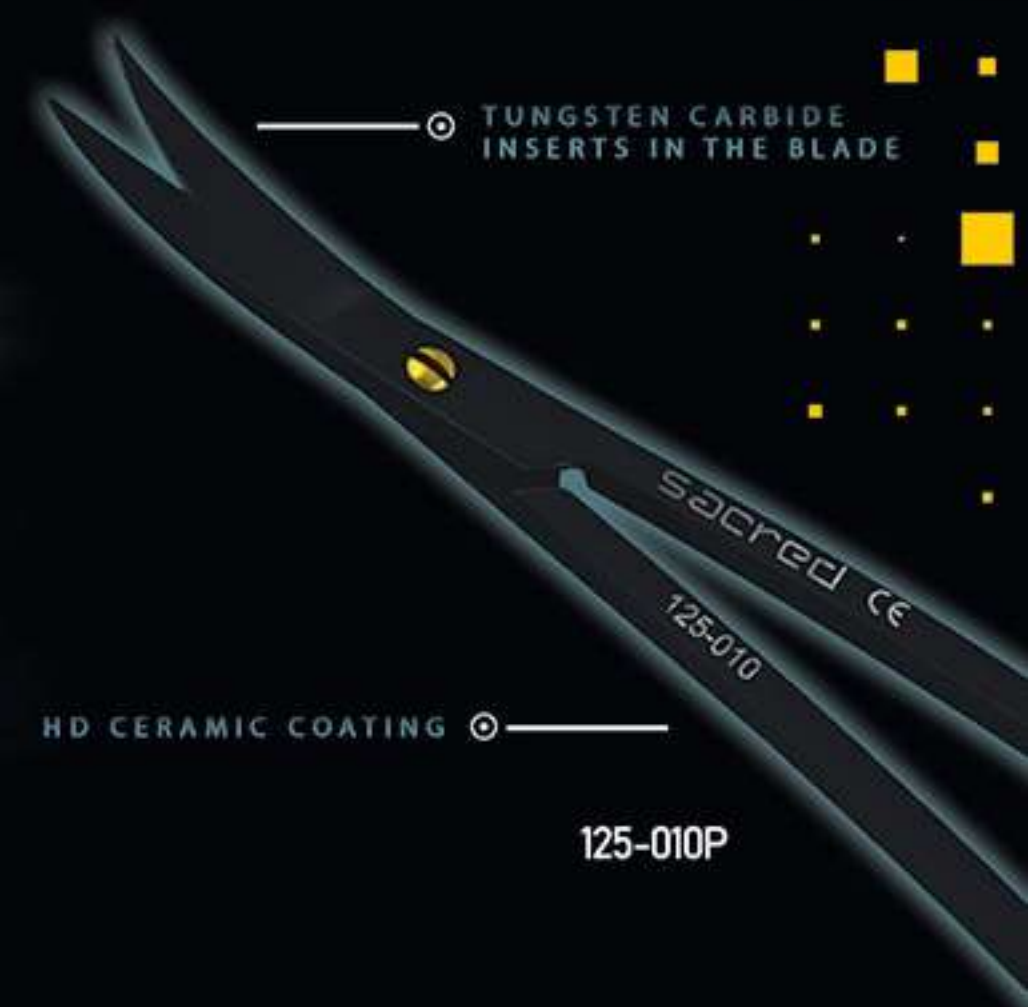

125-010P

COMPOSITE INSTRUMENTS (SET OF 6) - TITANIUM COATED



(SET OF 6) 025-600S

BONE GRAFTING INSTRUMENTS (SET OF 6) - TITANIUM COATED



(SET OF 6) 025-700S

www.sacred.dental

P.K. THOMAS INSTRUMENTS (SET OF 5) - STAINLESS STEEL



PKT-SET-SS



Autoclavable plastic instrument cassettes available in all above colours



PKT-SET-AL

* P.K Thomas Set is also available in aluminium, HD Black Coated.

HAND INSTRUMENTS BLACK TITANIUM COATED TIPS

MORE DURABILITY AND EFFICIENCY



025-501

025-502

025-503

025-504

025-505

025-506

OSTEOTOMES IMPLANTOLOGY SURGERY



700-134
2.8mm



700-132
3.3mm



700-131
3.8mm



700-133
4.2mm



700-135
5.0mm



ISO 9001
ISO 13485
FDA APPROVED



AISI 304 & 410
STAINLESS STEEL
MEDICAL GRADE



CORROSION RESISTANCE
UNDER TYPICAL ROOM
TEMPERATURE
ATMOSPHERIC
CONDITIONS.



UNIVERSAL STANDARD
DIMENSIONS & LENGTH
OF INSTRUMENTS




PRECISE & MECHANICALLY
ALTERED DENTAL TIPS
FOR EFFICIENT WORKING

sacred[®]
Symbol of quality

DENTAL CARE KIT PLAQUE REMOVER TOOLS



AMZ-SAC-01

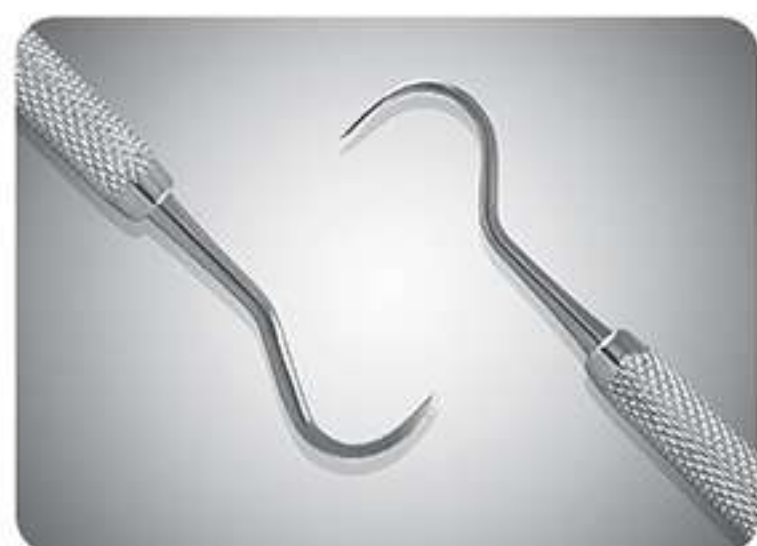

SICKLE SCALER

Use for tactile sensation to
detect tooth decay with
either side and to remove
food debris



DENTAL PROBE

Use to remove plaques
(which can't be removed
by toothbrush) on your
teeth



TARTAR SCRAPER

Use to remove tartar /
stains on your teeth
or your gum



DENTAL MOUTH MIRROR

To view the area of your
oral cavity that can't be
seen with a regular mirror

■ DENTAL CARE KIT PLAQUE REMOVER TOOLS

Our Plaque Remover kit is specifically designed and manufactured to remove dental calculus, plaque, tartar, food debris and stubborn stains, the mouth mirror No. #5 further helps you locate any tight invisible spots inside your mouth, hence guaranteeing overall hygiene satisfaction.



IMPLANTOLOGY KITS



EASY FAST S

This kit consist of:

- Torque ratchet 10 - 40 Ncm
- 4.0 mm square handle
- Hand wheels
- Triangle drill 1.8 mm
- Pilot 2.0 (white)
- Regular Drills
- Reamers

Regular Drills	Colour Indication	Actual Reamer Size	For Implant
2.5 mm	Purple	2.6 mm	3.0 mm
2.9 mm	Green	3.1 mm	3.5 mm
3.2 mm	Yellow	3.4 mm	3.8 mm
3.5 mm	Red	3.9 mm	4.3 mm
4.4 mm	Blue	4.6 mm	5.0 mm

DRILL GUIDE & IMPLANT POSITIONING KIT

The Implant Positioning Guide Kit is used for the placement of dental implants with the use of positioning guides and pins. The guides provided slide over a pilot drill which help center the drill between adjacent teeth. After the osteotomy is prepared, surgeons can preserve the space using a guide pin even when more than one implant is placed across the ridge making this kit an essential part of any implant practice.



THE OSSEODENSIFICATION BURS SURGICAL DRILLS BUR KIT

The osseodensification burs use external irrigation and are compatible with all available sized implants up to 6.0 mm leaving you with the flexibility to do more with less. The proprietary flute design allows the bur to work in two modes.

- 1). Cutting Mode (Clockwise Rotation)
- 2). Densifying Mode (Counterclockwise Rotation)

(At 800 - 1500 RPMs)



UNIVERSAL PROSTHETIC KIT IMPLANT DRIVER

Universal Prosthetic Driver Kit is used to remove and fasten implant prosthetics. The screwdrivers covers a vast range of implant screws with complete set of components for tightening clinical and prosthetic screws to the correct torque value.

Helps to easily deliver a hex driver to restore the position and tighten the screw including healing, cover and abutment screw without dropping off.




IK007

BUR DRILLS WITH STOPPER IMPLANT SURGICAL KIT

(Kit Includes 21 Pieces)

Stopper			
L3	L5	L7	L9
L10	L11	L12	L13

Drills Size (mm)			
2.0	3.0	4.0	5.0
2.3	3.3	4.3	5.3
2.5	3.5	4.5	5.5

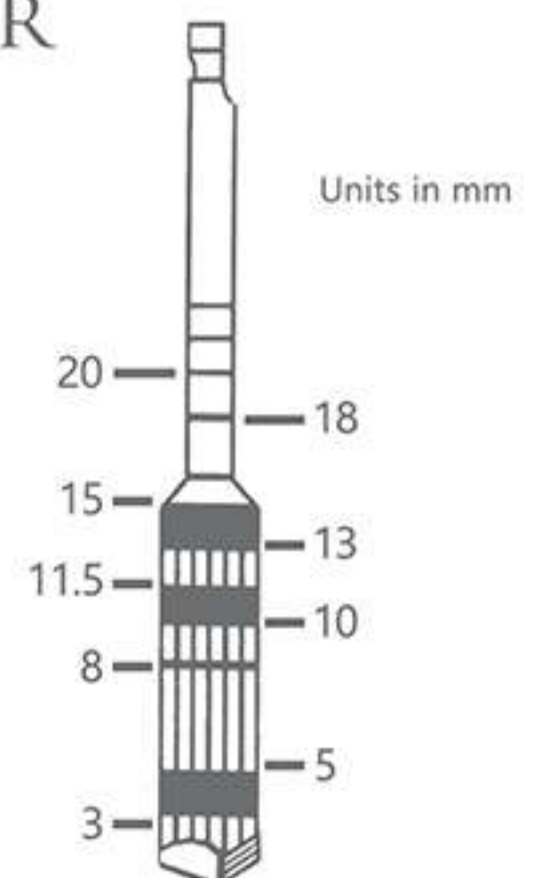

IK006



BONE EXPANDER KIT
IK001



SAW DISK KIT
IK004



SURGICAL DRILLS KIT
IK003



TISSUE PUNCH KIT
IK005

EXTRACTING FORCEPS



Available in different new patterns.

- (TiN) Titanium finish
- Impervious against rust
- Enhanced durability
- Offers aesthetic colours to meet your personal preferences.

MICRO SURGERY INSTRUMENTS



Castroviejo
Angled Point
14 cm



Castroviejo-Curved
19 cm
1915TC



Castroviejo-Straight
19 cm
1914TC



Castroviejo-Curved
NH551114.5 cm
NH561118cm



Needle holder
2015TC



Micro Forceps



Micro Surgical Scissor
Curved
18cm



Micro Surgical Scissor
Straight
18cm



Micro Surgical
Tissue Plier
18cm

SILICONE HANDLE INSTRUMENTS



025-410S

Available in different colours.



ALUMINIUM IMPRESSION TRAYS

sacred[®]
Symbol of quality



225-100SU



225-101MU



225-102LU



225-100SL



225-101ML



225-102LL

MOUTH MIRRORS



000-001-ES
Plane



000-017-ES
Front Surface



000-001-S
Plane



000-017-S
Front Surface



000-001-ES
Plane



000-017-ES
Front Surface



000-001-DS
Plane



000-017-DS
Front Surface

ORTHO INSTRUMENT



250-066

MINNESOTA SERRATED RETRACTOR




800-006 SER

Serrated end to fit the bone surface easily.



STERILIZATION CASSETTES

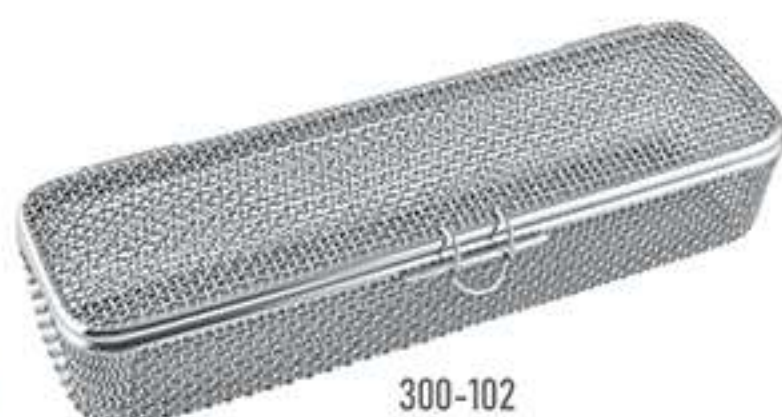

300-102



300-104



300-103



3D 300-105



3D 300-106

SCALPEL WITH DROP BUTTON



Scalpel with Drop Button plays a vital role by preventing contaminated injuries or infections. By pressing the button at the end of the scalpel handle the mechanism lifting the blade is triggered off. The blade is pushed forward so that it falls safely in a tray or cup.

After cleaning and sterilization, you can insert a new blade – as usual. The advantage lies in the fact users can remove the used blade in a safe and controllable way.



000-002 DB

ROOT ELEVATORS WITH (TiN) TITANIUM FINISH

Bein Root Elevator is used as a lever during the extraction of a tooth or retained roots. It is pushed in between the gums and the tooth exterior to loosen the tissue surrounding the tooth.



Bein
075-001BL

Cryer Elevator is used after coronal aspect of septum bone removal to luxate mandibular molar distal root in order to expose the distal root mesial surface. The triangular blade supports eradicating the broken roots in the tooth socket and when the adjacent socket is empty.



Cryer
075-044BL



Cryer
075-045BL

LABORATORY INSTRUMENTS



LA101

LA102


LA103

FAST ELASTIC POSITIONER

sacred
Symbol of quality

Fast Elastic Positioners are used for quick and efficient placement of elastomeric ligatures.

Engage the elastomeric by placing the tips in centre of the ring and pushing the lever back to expand. It can also be used to change the clear elastomeric ties when using ceramic brackets.

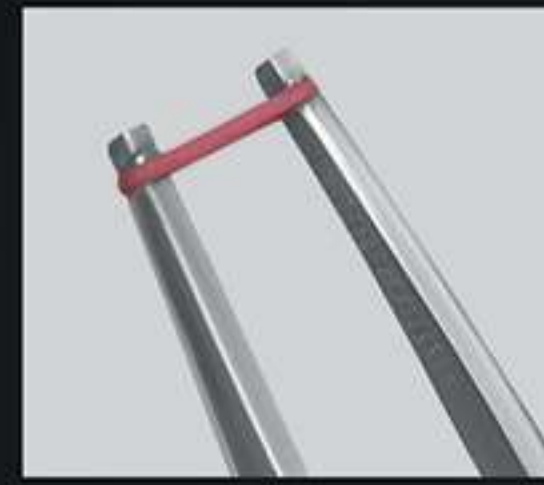

Upto 7.0 mm opening



Locking mechanism to expand elastomeric ligatures.



CAPSULE APPLICATOR BRIGHT

Capsule Applicator is used to extruding substances like sealants, bonding, composites and impression materials.



DENTAL RATCHETS



CORD PACK UNIVERSAL




IMPLANT SCALER



Note: Detailed catalogue is available on request.

www.sacred.dental

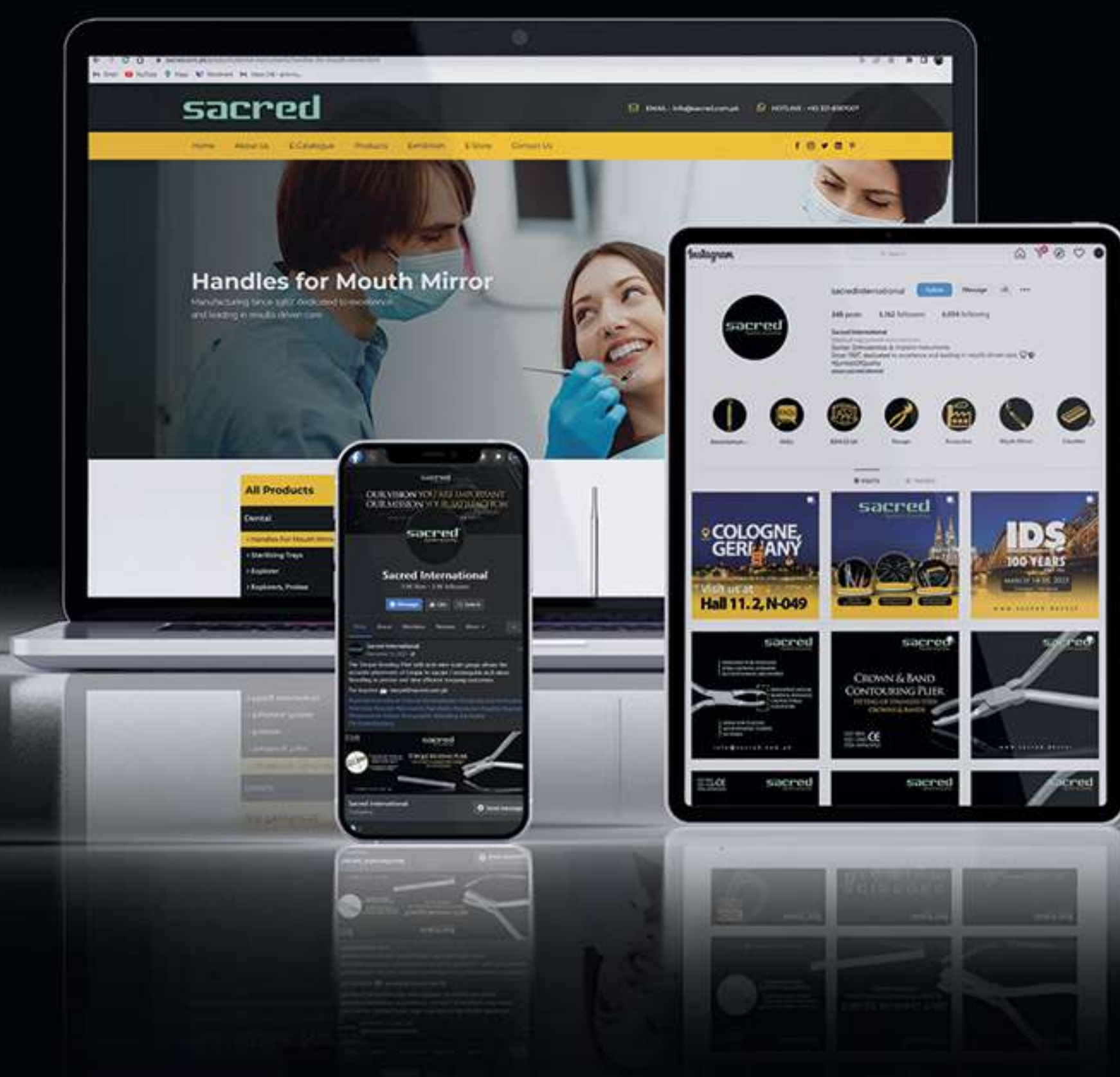

sacred[®]
Symbol of quality

📍 Harrar, Roras Road,
Sialkot-51310, Pakistan.

☎ +92 52 3557007
☎ +92 52 3257007

☎ +92 321 7227007

✉ info@sacred.com.pk

🌐 www.sacred.com.pk